


UNITY Aneuploidy NIPT: 22q11.2 Microdeletion Analysis

For all singleton and twin pregnancies



22q11.2 Deletion Syndrome is **Prevalent**

22q11.2 (also known as DiGeorge Syndrome) is more prevalent than commonly screened for chromosomal conditions.

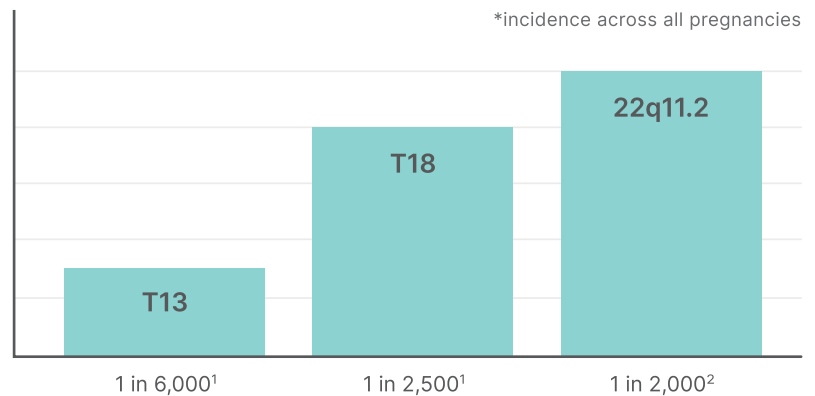

Early **Identification** is Difficult with Ultrasound

Ultrasound findings are often not unique to 22q11.2 and can make diagnosis difficult. 90% of 22q11.2 microdeletions are de-novo and prevalence does not increase with maternal age.³



Early Prenatal **Diagnosis** Can Make a Difference

Prenatal detection of 22q11.2 Deletion Syndrome is associated with less cardiac and noncardiac morbidity.⁴

Some babies can develop neonatal hypocalcemia that can cause seizures if not treated.⁵



Expertise in Fetal DNA
Quantification Allows for
Highly Accurate Detection
of 22q11.2 Microdeletion

Includes the full A-D region AND nested microdeletions.⁶

>95%
sensitivity

>99.9%
specificity

80%
PPV


Fetal DNA Quantification Overcomes the Technical Limitations with SNP Based Tests

QCT Counting Technology Allows for Extreme Precision and Sensitivity

We are experts at fetal DNA quantification, which allows for detecting small changes, even down to a single base pair.

QCT technology provides increased accuracy for smaller chromosomal changes, like 22q11.2.

This makes it possible to provide industry-leading sensitivity and PPV.

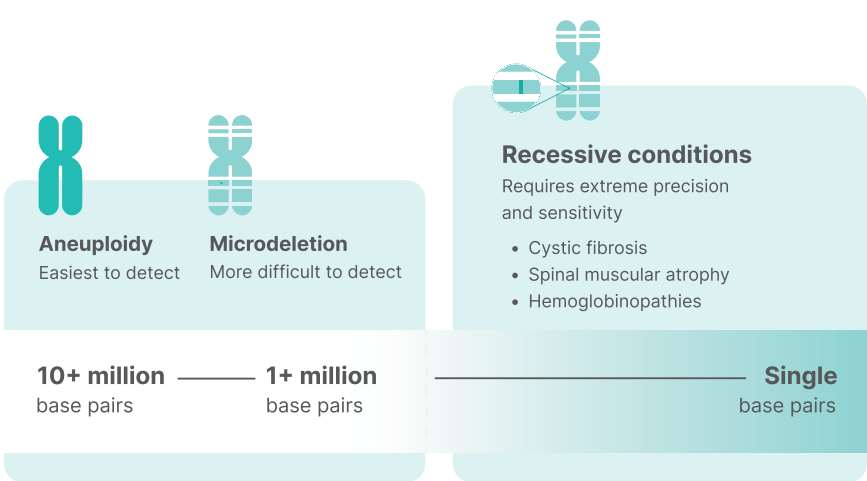

UNITY Aneuploidy Screen. Backed by clinical data.

New Publication: *Performance Characteristics of a Next Generation Sequencing-Based cfDNA Assay for Common Aneuploidies in a General Risk Population*⁷

	Trisomy 21	Trisomy 18	Trisomy 13	Combined Autosomes
Mean test characteristics n=114,707				
MATERNAL AGE 29 years old				
GESTATIONAL AGE 13.9 weeks				
FETAL FRACTION 9.30% (1.5–39%)				
TURNAROUND TIME 5 days				
Sensitivity	99.7%	99.5%	>99.9%	99.7%
Specificity	99.7%	>99.9%	>99.9%	99.9%
PPV	90.5%	97.6%	73.3%	90.8%
NPV	>99.9%	>99.9%	>99.9%	>99.9%

*1.5–39% represents the full distribution of fetal fraction

Also available with UNITY Aneuploidy Screen

Add-on at any point during the patient’s pregnancy

- + UNITY Fetal RhD NIPT
- + UNITY Fetal Antigen NIPT
- Twin zygosity determination*

*included with every twin order

Access the
peer-reviewed
publication here



650.460.2551 | support@unityscreen.com | unityscreen.com

2024 BillionToOne, Inc. All rights reserved. UNITY Fetal RhD is a trademark of BillionToOne, Inc. 1. Goel, Nitin et al. "Trisomy 13 and 18-Prevalence and mortality-A multi-registry population based analysis." American journal of medical genetics. Part A vol. 179,12 (2019): 2382-2392. doi:10.1002/ajmg.a.61365. 2. Fung, W., Butcher, N., Costain, G. et al. Practical guidelines for managing adults with 22q11.2 deletion syndrome. Genet Med 17, 599–609 (2015). https://doi.org/10.1038/gim.2014.175. 3. McDonald-McGinn DM, Hain HS, Emanuel BS, et al. 22q11.2 Deletion Syndrome. 1999 Sep 23 [Updated 2024 May 9]. In: Adam MP, Feldman J, Mirzaa GM, et al., editors. GeneReviews® [Internet]. Seattle (WA): University of Washington, Seattle; 1993-2024. Available from: https://www.ncbi.nlm.nih.gov/books/NBK1523/. 4. Freud LR, et al. Prenatal vs postnatal diagnosis of 22q11.2 deletion syndrome: cardiac and noncardiac outcomes through 1 year of age. Am J Obstet Gynecol. 2024 Mar;230(3):368.e1-368.e12. doi: 10.1016/j.ajog.2023.09.005. Epub 2023 Sep 16. PMID: 37717890. 5. Cheung, E. N., George, S. R., Andrade, D. M., Chow, E. W. C., Silversides, C. K., & Bassett, A. S. (2014). Neonatal hypocalcemia, neonatal seizures, and intellectual disability in 22q11.2 deletion syndrome. Genetics in Medicine, 16(1), 40–44. https://doi.org/10.1038/gim.2013.71. 6. UNITY Aneuploidy NIPT Validation. Data on File. March 2024. 7. Wynn J, Rego S, Fang J, Alford B, Carter R, Hoskovec J. Performance Characteristics of a Next Generation Sequencing-Based cfDNA Assay for Common Aneuploidies in a General Risk Population. Ann Gynecol Obstetr Res. 2024; 7(1): 1027. UN_FL_22Q_11-2_2407